

PREPARING THE INFUSION- USE ASEPTIC TECHNIQUE

Check the medicine label for:

- Your name
- The correct medication and dosage
- Expiration date

1. Wash your hands. If you have been taught to flush your IV catheter, do so at this time as instructed by your nurse.

Steps 2-5 explains how to prepare your infusion if you have a rubber cap on your IV catheter. If you have a Luer connector on your IV catheter, your nurse will provide other instructions.

2. Remove the cap at the end of the tubing and attach a needle or needleless device.

3. If you need to remove the air from the tubing, open the clamp. When the fluid has reached the end of the needle or needleless device and all air is removed from the tubing, close the clamp.

4. Wipe the rubber cap on your IV catheter with an alcohol pad. Allow the cap to dry completely.

5. Remove the needle or needleless device cover and insert the needle or connector into the center of the rubber cap on your IV catheter. Secure the connection with tape.

FOR C-SERIES HOMEPUMP USERS ONLY:

The lower portion of the tubing containing the flow restrictor should be worn in contact with your skin. Tape if necessary.

6. To start the infusion, open the clamp. If the tubing is kinked, massage the tubing to begin flow.

7. Place the pump in the carrying pouch, your pocket, or in a stationary place while getting your medicine.

END OF INFUSION

1. Your infusion is complete when the infusion time has passed and the balloon is no longer full. Close the clamp and remove the needle from the rubber portion of your IV catheter.

2. Dispose of the entire unit in an appropriate manner, as instructed by your nurse.

3. Flush your IV catheter as instructed by your nurse.

PRECAUTIONS

- Do not reuse.
- Do not use while bathing, showering, or swimming.

YOUR MEDICATION SCHEDULE:

Name of medicine:

Dosage:

Times to give medicine:

Length of medicine infusion:

IV catheter flush directions:

Medicine side effects to report:

If you have any questions about your therapy, contact:

PATIENT GUIDELINES

This guide will help you learn the simple steps in using the HOMEPUUMP Eclipse™ or the HOMEPUUMP “C” Series™ to give your medicine. The HOMEPUUMP is an easy to use IV delivery system.

Your nurse may give you more information about your medicine and IV therapy.

STORAGE

Your nurse will tell you if your medicine needs to be stored in a refrigerator, freezer, or at room temperature.

Your nurse will also tell you when to remove the HOMEPUUMP from the refrigerator or freezer before use. The chart in this guide will give you and your nurse this information.

Remove the HOMEPUUMP from your refrigerator or freezer _____ hours before use.

Do not use a microwave, oven, or bath-water to warm the fluid. Allow the HOMEPUUMP Eclipse to warm at room temperature. The HOMEPUUMP must be used at room temperature; the cold fluid will result in a longer infusion time.


ESTIMATED TIME REQUIRED TO REACH ROOM TEMPERATURE PRIOR TO USE

Model#	Description	Fill Volume	Flow Rate	Remove From Freezer	Remove From Refrig.
E050500	Homepump Eclipse	50 ml	50 ml/hr	12 hrs	6 hrs
E101000	Homepump Eclipse	100 ml	100 ml/hr	12 hrs	6 hrs
E100500	Homepump Eclipse	100 ml	50 ml/hr	12 hrs	6 hrs
E102000	Homepump Eclipse	100 ml	200 ml/hr	12 hrs	6 hrs
E251750	Homepump Eclipse	250 ml	175 ml/hr	18 hrs	10 hrs
E252500	Homepump Eclipse	250 ml	250 ml/hr	18 hrs	10 hrs
E401000	Homepump Eclipse	400 ml	100 ml/hr	24 hrs	12 hrs
E402000	Homepump Eclipse	400 ml	200 ml/hr	24 hrs	12 hrs
E502500	Homepump Eclipse	500 ml	250 ml/hr	30 hrs	18 hrs
C060020	Homepump C-Series	60 ml	2.0 ml/hr	12 hrs	6 hrs
C065005	Homepump C-Series	65 ml	0.5 ml/hr	12 hrs	6 hrs
C100005	Homepump C-Series	100 ml	0.5 ml/hr	12 hrs	6 hrs
C100020	Homepump C-Series	100 ml	2.0 ml/hr	12 hrs	6 hrs
C125050	Homepump C-Series	125 ml	5.0 ml/hr	12 hrs	6 hrs
C270010	Homepump C-Series	270 ml	1.0 ml/hr	24 hrs	12 hrs
C270020	Homepump C-Series	270 ml	2.0 ml/hr	24 hrs	12 hrs
C270050	Homepump C-Series	270 ml	5.0 ml/hr	24 hrs	12 hrs
C270100	Homepump C-Series	270 ml	10 ml/hr	24 hrs	12 hrs
C300060	Homepump C-Series	300 ml	6 ml/hr	24 hrs	12 hrs

CAUTION: Federal Law (USA) restricts this device for sale by or on the order of a Healthcare Professional.

HOMEPUUMP
ECLIPSE®

The Preferred Elastomeric
Infusion System

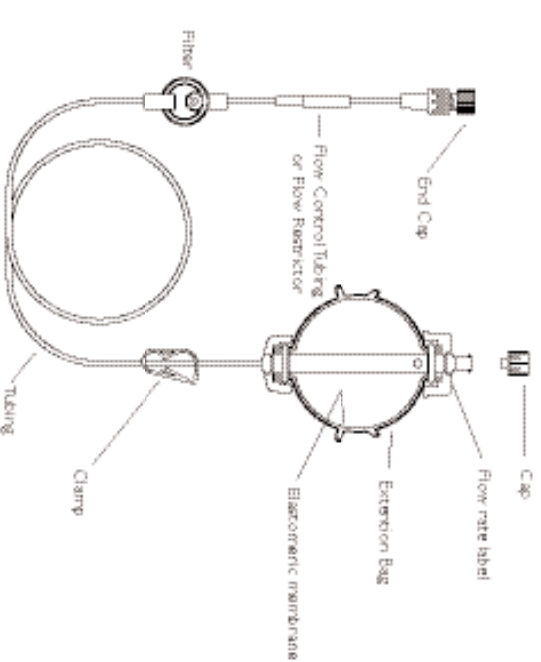

I-Flow Corporation
20202 Windrow Dr.
Lake Forest, CA 92630
Telephone: 800.448.3569
Clinical Hotline: 800.444.2728
Internet: www.iflo.com

HOMEPUUMP Eclipse® is a registered trademark of I-Flow Corporation
U.S. Pat. Nos. 5,080,652; 5,105,983; 5,284,481; and Foreign Pat. Pending
Printed in U.S.A